


Stevia glycosides

Stevia rebaudiana Bertoni is a plant used for its sweet taste. It includes glycosides, which are used as sweeteners in the food industry.

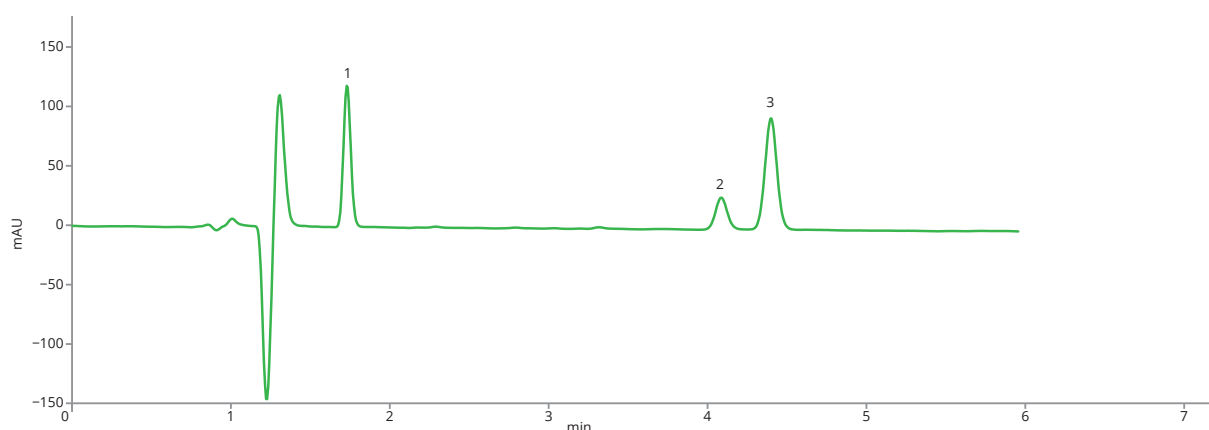

Standard on ASTRA® C18-HE

Column	ASTRA® C18-HE, 3 µm
Dimensions	150 mm × 4.6 mm
Part number	AST-5732-IK46
Mobile phase	Acetonitrile + 0.1% formic acid : water (32 % : 68 % v/v) Isocratic elution
Flow rate	1.6 mL/min
Temperature	65 °C
Detection	UV @210 nm
Analytes	1. Rebaudioside D 2. Rebaudioside A 3. Stevioside